One Component Self-Etching Light-Cured Dental Adhesive

**TOKUYAMA BOND FORCE II**

**TOKUYAMA BOND FORCE II** *Pen*

**CE** 1639 REF AD01

## Instructions for Use

Manufacturer

**Tokuyama Dental Corporation**

38-9, Taitou 1-chome, Taitou-ku, Tokyo 110-0016, Japan

Tel: +81-3-3835-7201

URL: <http://www.tokuyama-dental.com/>

EC Rep.

**Tokuyama Dental Italy S.r.l.**

Via Chizzalunga, 1

36066 Sandrigo, Vicenza, Italy

Tel: +39 0444 659650

Fax: +39 0444 750345

Importer

**Tokuyama Dental Deutschland GmbH**

Am Landwehrbach 5

48341 Altenberge Germany

Tel: +49 2505 938513

Fax: +49 2505 938515

<b>ENGLISH</b>	2	<b>ČESKY</b>	9	<b>DANSK</b>	16	<b>DEUTSCH</b>	24
<b>EESTI KEEL</b>	32	<b>FRANÇAIS</b>	40	<b>ITALIANO</b>	48	<b>LATVISKI</b>	55
<b>LIETUVIŲ K.</b>	63	<b>MAGYAR</b>	70	<b>NEDERLANDS</b>	77	<b>NORSK</b>	85
<b>POLSKI</b>	92	<b>SLOVENČINA</b>	100	<b>SUOMI</b>	107	<b>SVENSKA</b>	114

Read all information, precautions and notes before using.

## ■ PRODUCT DESCRIPTION AND GENERAL INFORMATION

1. TOKUYAMA BOND FORCE II is a one-component, one-coat application, self-etching, light-cured, dental adhesive system with excellent handling and adhesive properties and exceptional marginal integrity toward cut/uncut enamel and dentin when used in combination with light- or dual-cured composite materials.
2. A light-curing unit having the camphorquinone (CQ) wavelength range (peak: 470nm, spectrum: 400 to 500nm) can be used for curing TOKUYAMA BOND FORCE II.
3. TOKUYAMA BOND FORCE II contains Phosphoric acid monomer, Bisphenol A di(2-hydroxy propoxy) dimethacrylate (Bis-GMA), Triethylene glycol dimethacrylate, 2-Hydroxyethyl methacrylate (HEMA), Camphorquinone, alcohol and purified water. The pH level immediately after dispensing is approximately 2.8.
4. TOKUYAMA BOND FORCE II is available in three convenient dispensing systems, Bottle, Pen-Type and Unit Dose. Please check with your local distributor for availability.

## ■ INDICATIONS

Bonding of light- or dual-cured composite material to:

- cut/uncut enamel,
- cut/uncut dentin,
- fractured porcelain/composite repair.

## ■ CONTRAINDICATIONS

TOKUYAMA BOND FORCE II contains methacrylic monomers, organic solvents and acids. DO NOT use TOKUYAMA BOND FORCE II for patients allergic to or hypersensitive to methacrylic monomers, related monomers, organic solvents and acids.

## ■ PRECAUTIONS

- 1) DO NOT use TOKUYAMA BOND FORCE II for any purpose other than those listed in these instructions. Use TOKUYAMA BOND FORCE II only as directed herein.
- 2) TOKUYAMA BOND FORCE II is designed for sale and use by licensed dental care professionals only. It is not designed for sale nor is it suitable for use by non-dental care professionals.

- 3) DO NOT use TOKUYAMA BOND FORCE II if the safety seals are broken or appear to have been tampered with.
- 4) If TOKUYAMA BOND FORCE II causes an allergic reaction or oversensitivity, discontinue its use immediately.
- 5) Use examination gloves (plastic, vinyl or latex) at all times when handling TOKUYAMA BOND FORCE II to avoid the possibility of allergic reactions from methacrylic monomers. Note: Certain substances/materials may penetrate through examination gloves. If TOKUYAMA BOND FORCE II comes in contact with the examination gloves, remove and dispose of the gloves, and wash hands thoroughly with water as soon as possible.
- 6) Avoid contact of TOKUYAMA BOND FORCE II with eyes, mucosal membrane, skin and clothing.
  - If TOKUYAMA BOND FORCE II comes in contact with the eyes, thoroughly flush eyes with water and immediately contact an ophthalmologist.
  - If TOKUYAMA BOND FORCE II comes in contact with the mucosal membrane, wipe the affected area immediately, and thoroughly flush with water after the restoration is completed. Affected areas may whiten from protein coagulation, but such whitening should disappear within 24 hours. If such whitening does not disappear within 24 hours, immediately contact a physician, and the patient should be so advised.
  - If TOKUYAMA BOND FORCE II comes into contact with the skin or clothing, immediately saturate the area with alcohol-soaked cotton swab or gauze.
  - Instruct the patient to rinse his mouth immediately after treatment.
- 7) TOKUYAMA BOND FORCE II should not be ingested or aspirated. Ingestion or aspiration may cause serious injury.
- 8) To avoid the unintentional ingestion of TOKUYAMA BOND FORCE II, do not leave it unsupervised within the reach of patients and children.
- 9) DO NOT expose TOKUYAMA BOND FORCE II or the vapor to open flame because TOKUYAMA BOND FORCE II is flammable.
- 10) To avoid cross infection or lowering bond strength, DO NOT reuse the disposable applicator. Clean the dispensing well thoroughly with alcohol after each use.
- 11) DO NOT use TOKUYAMA BOND FORCE II directly to the cavity being in close proximity to the pulp. Pulp protection with glass ionomer lining or calcium hydroxide is recommended.
- 12) When using a light-curing unit, protective eye shields, glasses or goggles should be worn at all times.
- 13) Be aware that when a restoration is chipped from misaligned occlusion or bruxism (clenching,

grinding or tapping), the repaired restoration may chip again.

- 14) DO NOT mix TOKUYAMA BOND FORCE II with other brands of primers or adhesives.
- 15) TOKUYAMA BOND FORCE II should be brought to room temperature before using if it is stored under refrigeration.
- 16) If using the Pen dispensing system,
  - DO NOT push the button when the cap is closed.
  - DO NOT push the button when the nozzle is facing up.
  - DO NOT hold the Pen diagonally or horizontally while dispensing.
  - DO NOT use the Pen directly in the oral cavity
  - DO NOT dismantle the Pen.

## ■ PRECAUTIONS FOR MEDICAMENTS AND MATERIALS

- 1) Some materials and medicaments (hemostatic material) inhibit the adhesion of TOKUYAMA BOND FORCE II for an extended period even after careful meticulous cleansing with water. DO NOT USE products which contain:
  - eugenol,
  - diammine silver fluoride [molecular formula:  $\text{Ag}(\text{NH}_3)_2\text{F}$ ],
  - phenols such as parachlorophenol, guaiacol, phenol,
  - aluminum chloride,
  - ferric sulfate,
  - aluminum sulfate,
  - epinephrine.
- 2) TOKUYAMA BOND FORCE II does not adhere to tooth structure immediately after applications of materials which contain:
  - hydrogen peroxide (oxydol),
  - sodium hypochlorite

We advise avoiding the use of these materials on the day of cementation, however these materials can be used for root canal treatment as the inhibition caused by these materials usually evaporates within five days.

## ■ STORAGE

- 1) Store TOKUYAMA BOND FORCE II at a temperature between 0 to 25°C (32 to 77°F). Storage under inappropriate condition will accelerate the degradation of TOKUYAMA BOND FORCE II.

- 2) If not in use for a prolonged period of time, store in a cool dark place.
- 3) Keep away from heat, direct sunlight, sparks and open flames.
- 4) DO NOT use after the indicated expiration date on the bottle/package.

## ■ DISPOSAL

Unused TOKUYAMA BOND FORCE II is to be absorbed into an inert absorbent material such as gauze or cotton, and disposed of in accordance with local regulations.

## ■ CLINICAL PROCEDURES

### 1. Cleaning

Thoroughly clean the tooth surface with a rubber cup and a fluoride-free paste then rinse with water.

### 2. Isolation

The rubber dam is the preferred method of isolation.

### 3. Cavity preparation

Prepare the cavity and rinse with water. Add bevels to the enamel margins of anterior preparations (class III, IV, V), as well as chamfers to the margins of posterior preparations (class I, II). Bevels and chamfers assist in erasing demarcations between the cavity margins and the restoration, enhancing both esthetics and retention.

- If there is an adherent surface including uncut enamel, apply etching agent to the uncut enamel.
- In case of porcelain/composite repairs, roughen the surface with a bur or a diamond point to prepare the area for adhesion; apply phosphoric acid for cleaning; rinse thoroughly with water; air dry thoroughly and treat with a silane coupling reagent according to its manufacturer's instructions.

### 4. Drying

Dry the cavity by blotting or with an air syringe.

- DO NOT desiccate the vital tooth. Desiccation can lead to post-operative sensitivity.
- The substances listed below inhibit curing and should be removed from the tooth surface by thoroughly cleaning the tooth surface with alcohol, citric acid, or the application of phosphoric acid for 2 to 3 seconds before application:
  - 1) Oil mist from handpiece,
  - 2) Saliva, blood and exudates.

## **5. Pulp Protection**

Glass ionomer lining or calcium hydroxide should be applied if the cavity is in close proximity to the pulp. DO NOT USE EUGENOL BASED MATERIALS to protect the pulp as these materials will inhibit the curing of TOKUYAMA BOND FORCE II.

## **6. Dispensing**

### **6-1. Unit Dose**

Twist open the Unit Dose container and dip the disposable applicator.

- When opening the cap, make sure to verify the top and bottom of the container.
- To avoid cross infection, do not reuse the material on other patients.

### **6-2. Bottle**

Open the Bottle by removing the cap and dispense one or two drops of adhesive into the dispensing well. Immediately after dispensing, tightly replace the Bottle cap.


- Be sure to hold the Bottle of TOKUYAMA BOND FORCE II vertically when dispensing. Do not hold the Bottle diagonally or horizontally while dispensing; this could cause the adhesive to flow back into the nozzle and become contaminated.
- Wipe off excess adhesive from the tip of the nozzle prior to closing.

### **6-3. Pen**


Open the Pen by removing the cap when the cap is facing up and dispense one or two drops of adhesive into the dispensing well by holding the Pen vertically and pushing on button repeatedly until adhesive starts to flow (3 or 2 pushes, \* This depends on how the button is pushed and the position that the Pen is being held). Immediately after dispensing, tightly replace the Pen cap.

- Be sure to hold the Pen vertically when dispensing. Do not hold the Pen diagonally or horizontally while dispensing; this could cause the adhesive to flow back into the nozzle and become contaminated.
- Wipe off excess adhesive from the tip of the nozzle prior to closing.

## ■ How to open and close the Pen
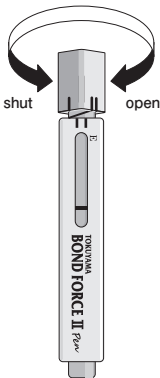
## ■ Indicator
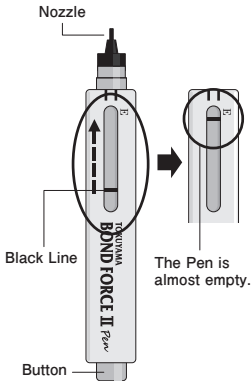
### Indicator:

The indicator allows users to monitor the remaining amount of adhesive. Black Line moves closer to the E as the amount of material decreases. When the line nears “E”, this indicates that the Pen is almost empty. This is not an exact reading but only an indication of the remaining material.

## 7. Application

Using the disposable applicator, apply TOKUYAMA BOND FORCE II on the surface to be bonded and wait 10 seconds or more.

- Ensure TOKUYAMA BOND FORCE II covers all surfaces where the composite resin will be applied.
- Protect the dispensed adhesive and the inserted applicator from ambient light before the application.

Unit Dose: keep the applicator in the Unit Dose container.

Bottle, Pen: use a light blocking plate.

- Complete the application within 5 minutes after dispensing because TOKUYAMA BOND FORCE II contains a volatile alcohol.
- Be careful not to allow excess TOKUYAMA BOND FORCE II to flow into the subgingival. If TOKUYAMA BOND FORCE II flows into the subgingival, thoroughly rinse and dry the area.
- If the adherent surface include ceramics, porcelain or composite, prepare the surface according to the appropriate procedures prior to the application of TOKUYAMA BOND FORCE II

## **8. Air Dry**

Using an oil-free air/water syringe, apply mild air continuously to the TOKUYAMA BOND FORCE II surface until the runny TOKUYAMA BOND FORCE II stays in the same position without any movement (usually for 5 seconds). Use a vacuum aspirator to prevent spatter of the adhesive.

- If accidental spattering occurs, it may cause the tissue to whiten or possible allergic reaction.
- If the adhesive pools on the cavity floor or the cavosurface angle and is too thick to air thin, blot the excess with a new disposable applicator before mild air application.
- If saliva, blood, or other fluids contaminate the applied adhesive, thoroughly rinse the cavity with water, dry and re-apply new adhesive.
- Do not rinse the applied adhesive with water except in unintentional contamination.

## **9. Light-cure**

Light-cure the surface for 10 seconds or more, keeping the curing light tip within a distance of 2 mm from the surface. If the cavity is too large or the light tip is too distant (e.g. MOD), divide the area into segments and light-cure each segment individually.

- Confirm that the light-curing unit has sufficient intensity ( $>300\text{mW}/\text{cm}^2$ ) before using. Note that using a cracked light guide will lower the intensity.

## **10. Light-cured Composite**

Restore with light-cured composite resin according to its manufacturer's instructions. The composite resin must be finished and polished.

- When dual-cured composite resins are placed into a cavity, the first increment must be light-cured using a layering technique.
- Do not use self-cured composite resins here because the Phosphoric acid monomer contained in TOKUYAMA BOND FORCE II may interfere with the curing of the self-cured resins, leading to their premature detachment.

The user and/or patient should report any serious incident that has occurred in relation to the device to the manufacturer and the competent authority of the Member State in which the user and/or patient is established.

**IMPORTANT NOTE:** The manufacturer is not responsible for damage or injury caused by improper use of this product. It is the personal responsibility of the user to ensure that this product is suitable for an appropriate application before use.

Specifications are subject to change without notice. When the product specifications change, the instructions and precautions may also change.

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