FOR CHAIRSIDE HARD DENTURE LININGS AND
EXTENSION OF DENTURE BORDERS

TOKUYAMA REBASE II

CE

Instructions for Use

Tokuyama Dental Corporation

38-9, Taitou 1-chome, Taitou-ku, Tokyo 110-0016, Japan

Tel: +81-3-3835-7201

URL: <https://tokuyama-dental.com>

ENGLISH	2	БЪЛГАРСКИ	10	ČESKY	19	DEUTSCH	27
DANSK	35	EESTI KEEL	43	ESPAÑOL	51	SUOMI	59
FRANÇAIS	67	ΕΛΛΗΝΙΚΑ	75	HRVATSKI	84	MAGYAR	92
ITALIANO	100	LIETUVIŲ K.	107	LATVISKI	115	NEDERLANDS	123
NORSK	131	POLSKI	139	PORTUGUÊS	147	ROMÂNĂ	155
SVENSKA	163	SLOVENČINA	171	SLOVENŠČINA	179	TÜRKÇE	187

ENGLISH

Read all information, precautions, notes and materials before using TOKUYAMA REBASE II.

■ Contents

Kit: TOKUYAMA REBASE II POWDER (80g)	× 1
TOKUYAMA REBASE II LIQUID (50mL)	× 1
ADHESIVE (15mL)	× 1
Measuring Cup (for TOKUYAMA REBASE II POWDER)	× 1
Dropper (for TOKUYAMA REBASE II LIQUID)	× 1
Rubber Cup (for mixing TOKUYAMA REBASE II POWDER and LIQUID)	× 1
Spatula (for spreading TOKUYAMA REBASE II on the denture)	× 1
Plastic Cup (for ADHESIVE)	× 2
Brush (for ADHESIVE)	× 2
TOKUSO RESIN HARDENER II	× 1
Spoon (for TOKUSO RESIN HARDENER II)	× 1

Refills: TOKUYAMA REBASE II POWDER (80g)
TOKUYAMA REBASE II LIQUID (50mL)
ADHESIVE (15mL)
TOKUSO RESIN HARDENER II

※Limited availability of FAST SET and/or NORMAL SET in certain areas. Contact your local dealer for availability.

■ Storage

1. Store the material at 0-25°C (32-77°F) at all times.
DO NOT use the material past its expiration date.
2. Replace caps tightly when finished.
3. If LIQUID is not to be used for more than a month, store the LIQUID with the original cap.
DO NOT store LIQUID with Dropper.
4. To avoid accidental ingestion or spillage, ALWAYS keep out of the reach of children and patients.

■ Disposal

1. To dispose of materials, mix the POWDER and LIQUID together, and allow it to harden. Be careful of the heat produced by the setting material.
2. Any remaining ADHESIVE should be absorbed in gauze or cotton, and disposed of in accordance with Federal and/or local regulations. Avoid any inhalation of the volatile substance.
3. Do not allow to enter drains, sewers or watercourses.
4. Must not be disposed of together with household garbage.
5. Do not rinse the POWDER bottle before disposal.
6. Follow local instructions for disposal.
7. Refer to the packaging waste disposal guide shown on the end of this IFU.

■ Contraindications

TOKUYAMA REBASE II contains methacrylic monomer. DO NOT use TOKUYAMA REBASE II for patients allergic to or hypersensitive to methacrylic monomers and related monomers or to any of the other ingredients.

■ Clinical Precautions

1. Use TOKUYAMA REBASE II only as instructed in these directions. DO NOT use TOKUYAMA REBASE II for any other application.
2. TOKUYAMA REBASE II is intended for use only by licensed dental care professionals.
3. Do not use this product if the safety seals are broken or if the contents have been tampered with.
4. Do not apply ADHESIVE to the patient's mouth or gums. Application of the ADHESIVE must be performed outside of the patient's mouth.
5. TOKUYAMA REBASE II LIQUID and ADHESIVE are flammable. Do not use near an open flame. Use in a well ventilated area.
6. Always clean the powder measuring cup, rubber cup, brushes, plastic cups and spatula with alcohol after each use.
7. Always determine whether your patient is allergic to, or is hypersensitive to acryl-based monomers. DO NOT use TOKUYAMA REBASE II on patients with a history of sensitivity to acrylic. Patients should be instructed to remove the prosthesis and immediately contact their dentist at the first sign of sensitivity.
 - Before using TOKUYAMA REBASE II, perform tissue conditioning for one week if there are signs of inflammation of the mucus membrane, or tissue irritaion from the old denture.

- Discontinue wearing the newly relined denture if there are any signs of hypersensitivity, redness, and ulcerations.
- 8. The use of examination gloves (plastic, vinyl, latex) is recommended at all times to avoid the possibility of allergic reactions from acryl-based monomers.
If ADHESIVE comes in contact with the examination gloves, immediately dispose of the examination gloves and wash hands thoroughly with water.
- After removal of the material from the patient's mouth, immediately have the patient rinse with water or mouthwash.
- 9. AVOID CONTACT WITH EYES. In case of contact with eyes, thoroughly flush eyes with water, and immediately contact a physician for assistance.
- 10. Avoid contact of uncured chemicals with skin. In case of contact with skin, clean affected area with rubbing alcohol, then thoroughly flush with water.
- 11. Immediately remove any excess lining material from the palatal area to avoid flow of material down the patient's throat. TOKUYAMA REBASE II is not radiopaque, and X-rays will not detect aspirated or ingested material.
- 12. Do not mix TOKUYAMA REBASE II with any materials not noted in these instructions.
- 13. To avoid breathing in the cured dust when polishing, always wear a mask recommended by Federal and/or local regulations.

■ Clinical Application

1. If the denture is a partial denture or has an undercut, residual root, etc, remove the denture from the mouth, trim it extraorally prior to final curing, and re-insert it for final curing. Otherwise, it will be difficult to remove the denture from the mouth.
 - FAST POWDER should only be used with FAST LIQUID, and NORMAL POWDER should only be used with NORMAL LIQUID.
 - Do not mix FAST POWDER with NORMAL LIQUID, or vice versa.
 - Measuring Cup for FAST POWDER is different from NORMAL SET.
 - Use corresponding Measuring Cup for each type.
2. Standard treatment time is based on temperatures of 23°C/73°F (Extraoral) and 37°C/98°F (Intraoral).
Treatment time is shortened by a rise in temperature caused by higher body temperature or elevated room temperature.
3. To minimize the formation of bubbles, mix TOKUYAMA REBASE II POWDER into the TOKUYAMA REBASE II LIQUID (White cap and White label).

4. Apply ADHESIVE (Blue cap and Blue label) only to denture surface. Do not mix ADHESIVE with TOKUYAMA REBASE II POWDER.

■ Patient Precautions

Tokuyama Dental recommends removing dentures at bedtime, in accordance with the following statement from the guidelines developed by the Global Task Force for Care of Full Dentures in August 2018:

“Denture wearers should not keep their dentures in the mouth overnight, unless there are specific reasons for doing so. This guideline is especially important for individuals at higher risk of developing stomatitis, as well as for frail or institutionalized older adults.”


■ Directions For Use

STANDARD POWDER-TO-LIQUID MIXING RATIO

One Graduation Powder

: Two Graduations Liquid


One Graduation
Powder

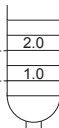Two Graduations
Liquid

Note: FAST and NORMAL come with different powder measuring cups.
Be sure to use each powder with its appropriate cup.

STANDARD TREATMENT TIME

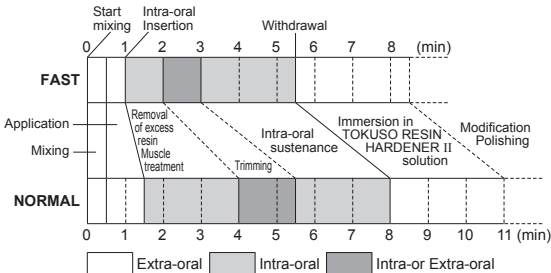*Treatment time varies according to room temperature and powder-liquid ratio.

- If trimming is done extraorally, always insert the denture back into the mouth just before final curing starts.
- Always keep the denture in the patient's mouth until final hardening is completed.
- Standard setting time is given at 23°C/73°F (extraoral) and 37°C/98°F (intraoral).

(1) Checking the adaptation

(2) Removal of the dentures surface

Carefully remove the old resin and expose the new resin of the dentures surface, especially the surface of the mucosal side, the marginal, buccal and labial areas where it is necessary to be relined with TOKUYAMA REBASE II.

Note:

- 1) Areas where the TOKUYAMA REBASE II should not contact, such as artificial teeth, the denture base or the interdental area, should first receive a coat of a separating agent.
- 2) As the relining surface of denture becomes larger, the resin paste shows a higher resistance in its flow. In this case, remove a large amount of the old resin or make a hole to let the resin flow out.

(3) Application of ADHESIVE

Apply ADHESIVE in a single layer with the brush to all areas of the surface to be relined. Leave the ADHESIVE for 20 seconds at room temperature. If the ADHESIVE liquid is still observed on

alveolar crest area of the denture, air blow to dry it completely.

Note:

- 1) Apply ADHESIVE to the buccal and labial side of the denture to be relined.
- 2) Apply ADHESIVE twice to the marginal area so that the relining material will not come off the denture.
- 3) ADHESIVE will remain effective for ten minutes after application.
- 4) Because ADHESIVE is extremely volatile, be sure to screw the cap tightly, immediately after use.

(4) Measuring and mixing

Use approximately one graduation of the TOKUYAMA REBASE II POWDER in the measuring cup for every two graduations of the TOKUYAMA REBASE II LIQUID in the dropper.

Note:

- 1) FAST and NORMAL come with different powder measuring cups. Be sure to use each powder with its appropriate cup.
- 2) The dropper for liquid can be used for either the FAST or NORMAL.
After measuring, place the powder into liquid (using the supplied rubber cup) and mix gently with the spatula for about 5-10 seconds. Mix gently to minimize bubbles.

Note:

- 1) Put the liquid in the cup first and then the powder to minimize the formation of bubbles.
- 2) Mixing too vigorously or for too long of a period will facilitate bubble formation, producing an uneven consistency. This in turn will create a rough surface and result in incomplete curing.
- 3) If necessary or so desired, the viscosity of the resin can be increased by increasing the quantity of the powder used.

(5) Application of the resin

With the spatula, apply the resin to the surface of the denture that will come in contact with the oral mucosa.

Note:

- 1) Apply resin within 20 to 60 seconds (at normal temperatures) after mixing for FAST (for NORMAL, within 30 to 90 seconds).
- 2) Apply the resin paste to all surfaces being relined, especially to the marginal, buccal and labial areas in the denture.
- 3) Be careful not to apply too much resin to the surface. Applying too much resin to the surface will result in the resin flowing down the patient's throat upon insertion of the denture.
- 4) Viscosity and hardening are temperature dependent. The warmer the room temperature, the

faster the set. The colder the room temperature, the slower the set.
If too much time is taken to mix the powder and liquid, the available working time will be shortened.

(6) Insertion of the denture

Insert the denture properly into the patient's mouth and have them bite down correctly at the centric position.

Note:

- 1) Insert denture within 60 to 90 seconds after the mixing has begun for FAST (for NORMAL, within 90 to 120 seconds).
- 2) Be extremely careful to avoid occlusal malposition.
- 3) Make sure there is no extra resin flowing down the patient's throat.
- 4) In case excess resin flows into the patient's throat, remove it with a fingertip or an appropriate instrument.

(7) Muscle Treatment and trimming

After following the above steps, and when the resin begins to harden, perform the necessary treatment of the muscles in the mouth. When the resin becomes harder than paste, remove the denture from the patient's mouth, and trim if necessary.

Note:

1) Timing

Be sure to complete all of the muscle treatment and trimming within 3 minutes from the time mixing started for FAST (within 5 1/2 minutes for NORMAL).

2) In case of a Full Denture

There is no need to remove the denture from the mouth, unless there are undercuts. Remove the excess portion of the resin intraorally, and keep the denture in the patient's mouth until hardening is complete.

3) In case of Undercut & Partial Denture

In the situation where patient has an undercut, partial denture, residual root or other problems, remove the denture from the mouth and trim it extraorally prior to final curing. Otherwise it will be difficult to remove the denture from the mouth.

(8) Re-insertion and oral maintenance of the denture (final curing)

Re-insert the denture and keep it in the patient's mouth until final hardening is complete.

Note:

- 1) Keep the denture fixed in the patient's mouth for 3 to 5 1/2 minutes after mixing has commenced for FAST. For NORMAL keep the denture fixed in the patient's mouth for 5 1/2 to

8 minutes. The resin will cure at this time. Do not remove the denture from the mouth during the period, since the presence of oxygen prevents the resin from fully curing.

- 2) Incomplete cure as a result of removing the denture too early will create undesired problems such as an inadequate surface hardness or a poor adaptation.

Using TOKUSO RESIN HARDENER II

1. For additional hardness of cured denture reline, TOKUSO RESIN HARDENER II maybe used.
2. TOKUSO RESIN HARDENER II improves the surface hardness and final polish of the denture base.
3. Prepare 40-60°C (104-140°F) water in a container.
Be sure to pour enough water to completely immerse the denture.
4. Dissolve 1 scoop of TOKUSO RESIN HARDENER II per 200mL of water.
5. Once the material is completely dissolved, immerse the denture for 3 minutes.
6. Remove the denture, rinse, and polish.

(9) Modifying and polishing

After final hardening, use a carbide bur, silicone point, finishing wheel, etc. to modify and polish the denture near the local exhaust ventilation.

(10) Checking the adaptation

(11) Completion

The user and/or patient should report any serious incident that has occurred in relation to the device to the manufacturer and the competent authority of the Member State in which the user and/or patient is established.

IMPORTANT NOTE: The manufacturer is not responsible for damage or injury caused by improper use of this product. It is the personal responsibility of the dental professional to ensure the product is suitable for application before use.

Specifications are subject to change without notice. When the product specification changes, the instructions and precautions may change also.

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Packaging waste disposal guide

(According to European Commission Decision 97/129/EC)

ENGLISH

Packaging waste disposal guide

Bottle (cap and body) for POWDER	Dispose of as industrial waste.
Bottle (cap and body) for LIQUID	Dispose of as industrial waste.
Bottle (cap and body) for ADHESIVE	Dispose of as industrial waste.
Bottle (cap and body) for RESIN HARDENER II	Dispose of as industrial waste.
Packaging bag for dropper etc.	LDPE4 (Low density polyethylene):Recycling. Follow local instructions for disposal.
Packaging bag for brush	PP5 (Polypropylene):Recycling. Follow local instructions for disposal.
Partition	PAP20 (Corrugated fibreboard):Recycling. Follow local instructions for disposal.
Paperbox	PAP21 (Non-corrugated fibreboard):Recycling. Follow local instructions for disposal.
IFU and guide	PAP22 (Paper):Recycling. Follow local instructions for disposal.

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